

CELLNERGY WELLNESS



INSTRUCTION MANUAL

[ENGLISH](#)

CLICK BELOW FOR OTHER LANGUAGES

[DEUTSCH](#) | [MAGYAR](#) | [POLSKI](#) | [ROMÂNĂ](#)



PRODUCT SPECIFICATION

Voltage: 110-127V for North America

220-240V for EU

Frequency: 50/60Hz

Power: 900W

Sound: 80dB

SAFETY INSTRUCTIONS

Please read the instruction manual carefully before using.


- Keep device away from water.
- Please supervise children under 12 years old or anyone with physical or mental impairments who are using, cleaning, and maintaining this device. Do not leave children unattended with device or allow children to play with device.
- Do not use device near baths, showers, sinks, or near any liquids.
- Please be sure to unplug the device after turning it off after each use, especially in bathrooms.

WARNING: ELECTROCUTION HAZARD

Do not submerge this device in water or other liquids. If the device comes into contact with water, immediately disconnect it from any power source and do not touch it. Contact a qualified electrician for assistance. Failure to follow these instructions may result in electrocution or other serious injury. Keep this device away from children and water sources and use only in a dry environment. Do not use this device if it has been damaged by water or other liquids.



DANGER: KEEP AWAY FROM WATER!

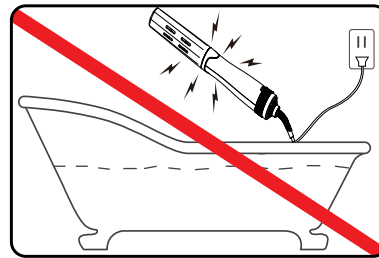

CAUTION

1. Read instruction manual carefully before use.
2. Do not allow children to play with device.
3. Please avoid knocking or dropping device.
4. Do not use the device outdoors, it should only be used in households.
5. Remove contact lenses before using device on face.
6. Do not plug in device if setting is not switched to OFF beforehand.
7. Do not blow device for more than 3 seconds on skin without movement.
8. Do not block the airflow of the device.
9. Do not use the device for more than 15 minutes daily.
10. Turn the setting to OFF before unplugging the device after each use.
11. Do not immerse the cord, plugs, or device in any liquids.
12. Do not place or store device where it can be fallen or drop into a tub or sink.
13. Do not reach for an appliance that has fallen into water. Unplug it immediately.
14. Do not use while bathing or in the shower.
15. Do not use the device if the supply cord is damaged.
16. If the supply cord arrived damaged, please contact our customer service team.
17. Do not leave the device unattended.
18. Unplug this appliance before cleaning.

SAFETY GUIDE

- Keep the product out of reach of children.
- Close supervision is necessary when the appliance is used by, on, or near individuals with certain disabilities or special needs.
- Do not leave appliance unattended when plugged in.
- Do not disassemble or modify the product in any way.
- Do not use the product near water or in wet conditions.
- Do not place the product near an open flame or any sources of heat.
- Do not expose the product to extreme temperatures or direct sunlight.
- Make sure the voltage of the product corresponds to the voltage of your power outlet before plugging in the device.
- Only use the product with the supplied power adapter.
- Do not use the product if the power adapter or cable is damaged, if it has been dropped or damaged, or if it has been submerged in water.
- Disconnect the device by gently pulling the plug of the power cord, not by pulling on the power cord directly.
- Do not use the product while driving or operating heavy machinery.
- Do not use the product if you have a pacemaker or any other medical device that could be affected by electromagnetic fields.
- Disconnect the device from the power socket after each use.
- Do not repair or make any adjustments to the device.
- Keep cord away from heated surfaces. Do not wrap the cord around the appliance.
- Do not use while sleeping.
- Do not operate when aerosol (spray) products are being used, or where oxygen is being administered.
- Do not use extension cord.
- Always store device in moisture-free area.
- Do not store device in temperatures exceeded 140°F (60°C).
- Do not unplug or plug this device when your hands are wet.

USER GUIDELINES

1. Do not use this product if you have hemorrhagic disease, a serious infection, hypertension, or open wounds.
2. If you have heart disease or any other medical condition, please check with your physician before using this product.
3. Do not use this product if you are undergoing surgical operations or pregnancy.

Disclaimer: This product is used for self-education and research purposes to support homeostasis. It is not intended to treat disease or to prevent impairment of human health. When used according to directions, it does not present any potential risk of illness or injury.

INSTRUCTIONS

1. Familiarize yourself with the settings:

Setting 0: Off

Setting 1: High speed, low heat

Setting 2: Low speed, high heat

Setting 3: Medium speed, high heat



We recommend Setting 1 for facial treatments, Setting 2 for ribcage and spine for stem cell activation, and Setting 3 for targeted areas.

2. Make sure the setting is OFF before plugging the device to the power outlet.
3. Plug the device into a power outlet and then use the dial at the bottom to adjust the setting.
4. Hold the device 8 to 12 inches (20 to 30 centimeter) away from the body. You can either blow the hot air in sweeping back and forth motions on a large area over bare skin or over clothing. Spend 10 minutes on each treated area on the body.
5. When you have finished using the device, turn it off before unplugging it from the power outlet.

MAINTENANCE

1. Always unplug and allow device to cool before cleaning.
2. Periodically inspect the air filter screen on the back of the brush and remove any debris that might have accumulated. Any build up may reduce air flow and damage the brush.
3. Never use thinner, gasoline alcohols, or any flammable chemicals to clean the product.
4. Clean the device with a soft dry cloth to remove any residue.
5. Do not wrap the cord around the device, this might cause the cord to wear out.

WARRANTY


One year warranty (from date of purchase): If the device is defective during the warranty period, please contact our customer service team.

Warranty will be void if:

1. Warranty period is exceeded.
2. The device is damaged due to not following the instruction manual.
3. The device is damaged due to mishandling or repair by an unauthorized party.

STEPS OF USING CELLNERGY WELLNESS

- Drink water before and after using the device.
- Hold the device 8 to 12 inches (20 to 30 centimeter) away from the body.
- Blow the hot air in sweeping back and forth or circular motions over bare skin or over clothing. Spend 10 minutes on each treated area on the body.

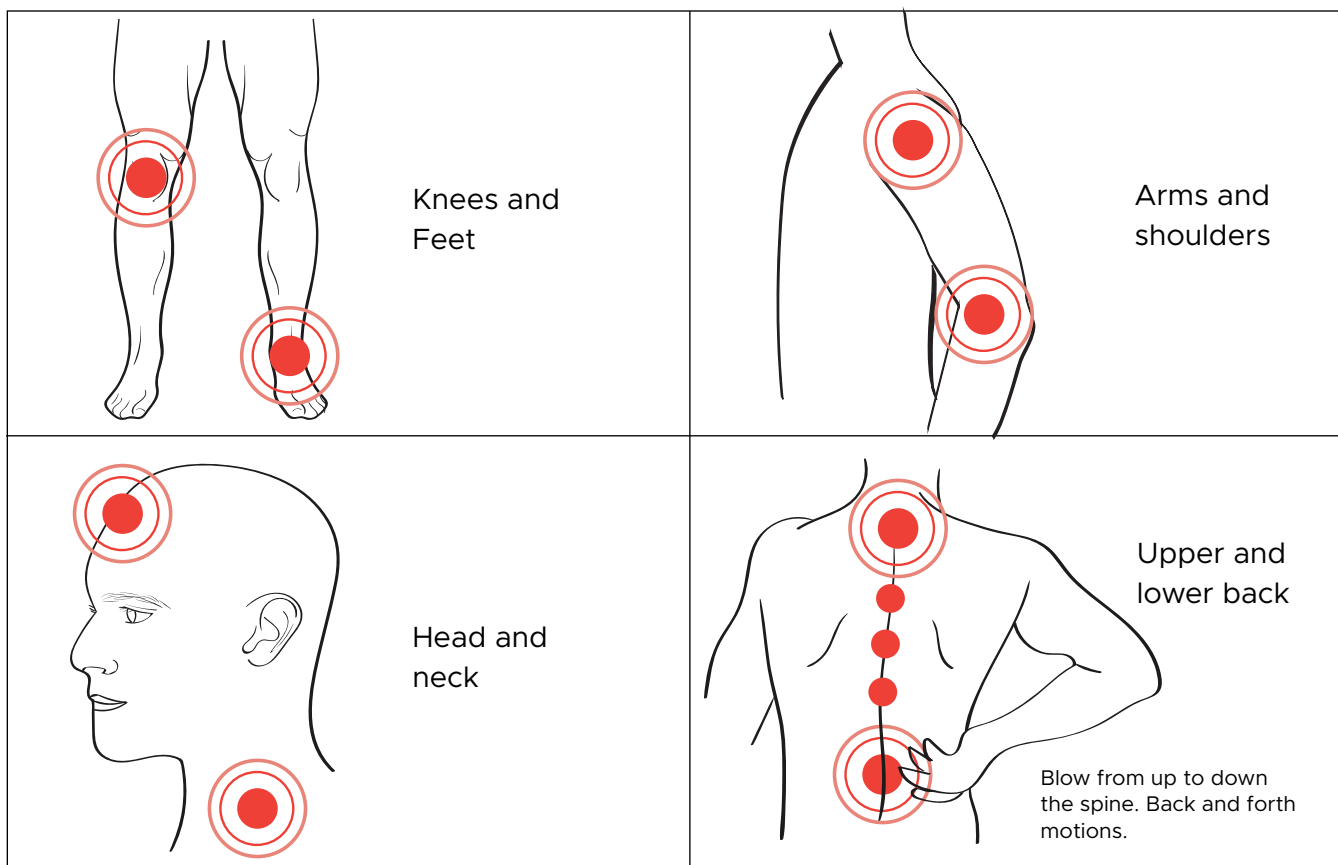

We recommend Setting 1 for facial treatments, Setting 2 for ribcage and spine for stem cell activation, and Setting 3 for targeted areas.

For best results, we recommend using Cellnergy Wellness daily with LAMININE and OMEGA⁺⁺⁺ supplements.