

CELLNERGY FOOTCHARGE



INSTRUCTION MANUAL

[ENGLISH](#)

CLICK BELOW FOR OTHER LANGUAGES

[DEUTSCH](#) | [MAGYAR](#) | [POLSKI](#) | [ROMÂNĂ](#) | [CZECH](#)



PRODUCT SPECIFICATION

Model: CFC5

Dimensions: 30 x 30 x 3 cm / 11.81 x 11.81 x 1.18 inches

Power: 5W

Typical Voltage in Use: 5V (USB 2.0 standard)

Typical Max Current: 0.5A

Typical Power Output: 2.5 watts

Temperature: Mild Warming (40-45 °C/104-113 °F)

Terahertz Chip: Embedded

Contains Lithium-Ion Battery

WHAT'S INSIDE



Cellnergy FootCharge
Device




Remote Control
(CR2032 Battery Installed)



USB
Charging Cable

PRODUCT CONFIGURATION

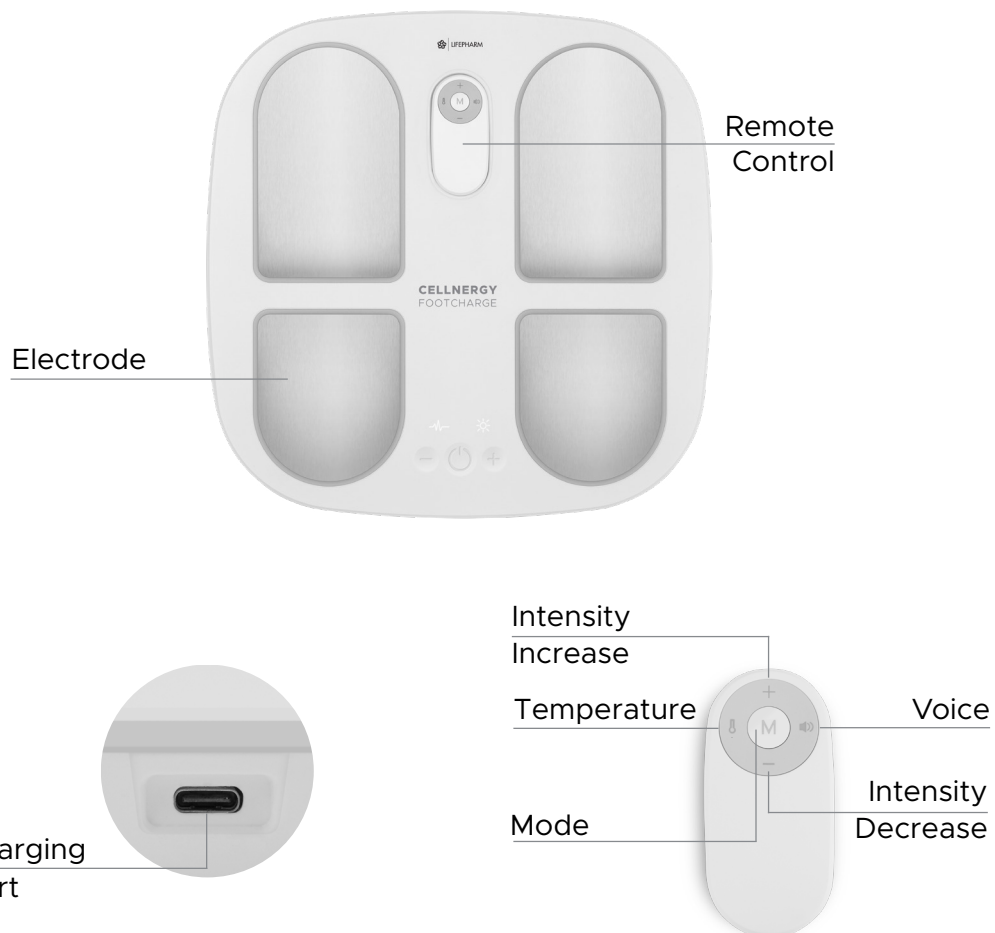

SAFETY INSTRUCTIONS

Please read the instruction manual carefully before using.

1. Use a damp towel to moisten your feet before use to improve conductivity.
2. This product is not waterproof. Do not immerse in water or use in damp areas.
3. Do not move the device during use. Pause or turn off the device before charging.
4. Do not allow children or individuals with disabilities to operate the device without supervision.
5. If discomfort occurs, stop using immediately and consult a physician.
6. Remove the battery if the product will not be used for a long time.
7. Keep metallic objects away from the electrodes during use.
8. Do not apply to the same area for extended periods.
9. Those sensitive or allergic to heat should consult a doctor before use.
10. Do not use if you have impaired skin sensation.
11. Avoid high-temperature or damp environments.
12. Do not use with medical implants such as pacemakers or ECG monitors.
13. Do not use simultaneously with high-frequency surgical devices.
14. Do not disassemble, repair, or modify the device yourself.
15. Do not use for commercial treatment purposes.

WARNING

- Keep out of reach of children.
- Do not incinerate.
- Do not expose to high temperatures.
- Do not immerse device in water.
- This is a personal health device, not a medical product.

WARNING: ELECTROCUTION HAZARD

Do not submerge this device in water or other liquids. If the device comes into contact with water, immediately disconnect it from any power source and do not touch it. Contact a qualified electrician for assistance. Failure to follow these instructions may result in electrocution or other serious injury. Keep this device away from children and water sources and use only in a dry environment. Do not use this device if it has been damaged by water or other liquids.

DANGER: KEEP AWAY FROM WATER!



BATTERY SAFETY

- Use only CR2032 batteries.
- Remove battery if not in use for extended periods.
- Do not use any other batteries other than the specified type.
- Do not use the battery when the temperature exceeds 45 °C/113 °F
- Dispose of batteries properly according to local laws.


USER GUIDELINES

1. This product is intended for people over the age of 16 only.
2. If you have any medical condition or concern, please consult with your physician before using this product.
3. Do not use if you have a pacemaker or implanted medical device.
4. Do not use if you are pregnant or menstruating.
5. Do not use if you have malignant tumors.
6. Do not use if you have foot wounds, infections, or circulatory disorders.
7. If you feel pain, dizziness, or chest discomfort during use, stop using the device immediately.
8. Do not use when you are under the influence of alcohol or sedatives.

Disclaimer: If you have any medical concerns or conditions, please consult with a physician before use. This product is intended for educational and research purposes to support general wellness. It is not designed to diagnose, treat, cure, or prevent any disease. When used as directed, it poses no risk of illness or injury.

INSTRUCTIONS

1. **SIT DOWN AND TAKE OFF YOUR SHOES AND SOCKS. PLEASE REMAIN SEATED THE ENTIRE TIME WHEN USING THE DEVICE, DO NOT STAND OR MOVE YOUR FEET.**
2. Place both feet on the electrodes at the same time. Ensure the skin of your feet is in direct contact with the electrodes.
3. Long press the “⏻” button on the device for 3 seconds. Follow the voice prompt. The device starts in full-function mode. Long press again to power off.
4. Press “M” to switch between modes. Press and hold the “M” button for 3 seconds to turn off.
5. Press “🌡️” to adjust heating temperature.
6. Press “🔊” to toggle voice guidance on or off. If you choose to leave the voice function on, please follow the voice prompt for guidance.
7. Press “+” or “-” to adjust electric pulse intensity.
8. The device automatically turns off after 15 minutes. It is recommended to use the device 15–30 minutes per session, up to twice daily. Do not exceed 60 minutes per day.
9. Turn off the device if you are done with use, or the device will automatically shut off after 15 minutes. Wipe clean with a soft dry cloth.

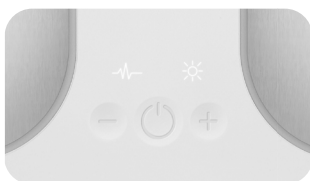

Long press “⏻” button for 3 seconds to turn the device on, and start to operate according to the voice prompt.
The default is full-function mode, flip to switch other modes.
Long press for 3 seconds again to turn it off.



Press “M” to switch modes
Press and hold the “⏻” button for 3 seconds to turn off.
Press “🌡️” to adjust temperature
Press “🔊” to turn on/off voice
Press “+” “-” to adjust pulse intensity

FOOTCHARGE MODES

Choose the mode that matches your body's needs:

MODE 1: RELAX

Supports deeper sleep and relieves nerve fatigue*

MODE 2: KNEAD

Boosts circulation and eases joint/muscle tension*

MODE 3: TAP

Improves balance, coordination, and reaction speed*

MODE 4: SQUEEZE

Enhances brain-body communication for clarity and focus*

MODE 5: PULSE

Eases numbness, tingling, and neuropathy-related symptoms*

MODE 6: STRENGTH

Supports muscle tone and helps prevent atrophy*

*These statements are intended to support general wellness and are not intended to diagnose, treat, cure, or prevent any disease.

TIPS

1. Please follow the voice guidance after turning on the device.
2. Always store the device in a dry and ventilated place after each use.
3. Do not put heavy objects on the device during use.
4. Please place both feet on the electrodes at the same time.
5. The device automatically turns off after 15 minutes.

WARRANTY

One-year warranty (from date of purchase) for product replacement. If the device is defective during the warranty period, please contact our Customer Service team through our website www.lifepharm.com or by emailing customerservice@lifepharm.com for a new or refurbished replacement.

Warranty will be void if:

1. Warranty period is exceeded.
2. The device is damaged due to not following the instruction manual.
3. The device is damaged due to mishandling or repair by an unauthorized party.

MAINTENANCE

1. Please use a dry soft cloth or towel to gently wipe the device after each use.
2. Please clean the electrodes after use with a damp soft cloth or towel.
3. Please use a small amount of rubbing alcohol with a concentration of 75% to disinfect the device with a soft cloth occasionally.

STORAGE

1. Please keep out of the reach of children.
2. Do not store in direct sunlight, high temperature, or damp places.
3. Store in a dry and ventilated place.
4. Do not store heavy items on top of the device.

FAQS

Q: Why do I feel tingling when using the electric pulse function

A: The tingling is from the CellWave BioActivation. If the intensity of the tingling is too high, please reduce the intensity by pressing the “-“ button.

Q: Why does the device feel too hot during use?

A: The heating function may take time to adapt. Reduce duration, lower the temperature, or stop use if discomfort continues.

Q: How long does it take to charge and how long does it last?

A: Charging takes about 3 hours. One full charge lasts approx. 3 uses (15 min twice a day).

Q: Why is there no noticeable increase in stimulation when intensity is increased?

A: Ensure both feet are damp and are in full contact with both electrodes.